

# Using the Calibration Block and the spectrum\_w\_cal.grc Program

## Instructions for Teachers and non-Gnuradio Experts

### Introduction

The horn telescope detects all signals incident on the antenna. We are primarily interested in the signal that is coming from the astronomical object of interest, such as the galaxy. So we need to subtract out any signals coming from any local terrestrial sources and from the telescope system. The **spectrum\_w\_cal.grc** Gnuradio program enables the user to display calibrated spectra with the unwanted signals subtracted out, while doing a run. This is accomplished through the **system\_calibration block**. This program also enables the user to “capture” the currently displayed spectrum to a .csv file (same as a .txt file), which can easily be opened in a spreadsheet.

### The system\_calibration Block

All signals detected by the horn antenna get amplified by the LNA. So the signal detected can be described by the following relation:

$$P_{\text{detected}}(f) = G(f) \times [T_{\text{object}}(f) + T_{\text{system}}(f)]$$

$P_{\text{detected}}(f)$  = total output power of the system at each frequency

$f$  = frequency;

$G(f)$  = gain at each frequency;

$T_{\text{object}}(f)$  = signal (in Kelvin) of the astronomical object of interest (what we want)


$T_{\text{system}}(f)$  = extraneous signals from the system and the ambient terrestrial environment

By measuring the signal from objects at known temperatures, the so-called “hot” and “cold” sources, the gain  $G(f)$  and the system temperature,  $T_{\text{system}}(f)$ , can be determined. Then on any subsequent measurements the signal object,  $T_{\text{object}}(f)$ , can be extracted.

The **system\_calibration** block does the calculations and outputs the  $T_{\text{object}}$  signal. This is done after the user completes the hot and cold calibrations.

## Features of the *spectrum\_w\_cal.grc* Gnuradio Program

- Data is calibrated and displayed on the fly without the need for any post-run manipulations of the spectrum and calibration data.
- During a run various display options are available: spectrum not calibrated, calibrated spectrum, hot calibration, and cold calibration.
- The scale of the signal axis can be adjusted during the run.
- The current spectrum being displayed on the main graph can be saved in a .csv text file for easy access.
- Writing to an hdf5 file can be stopped and continued during the run.

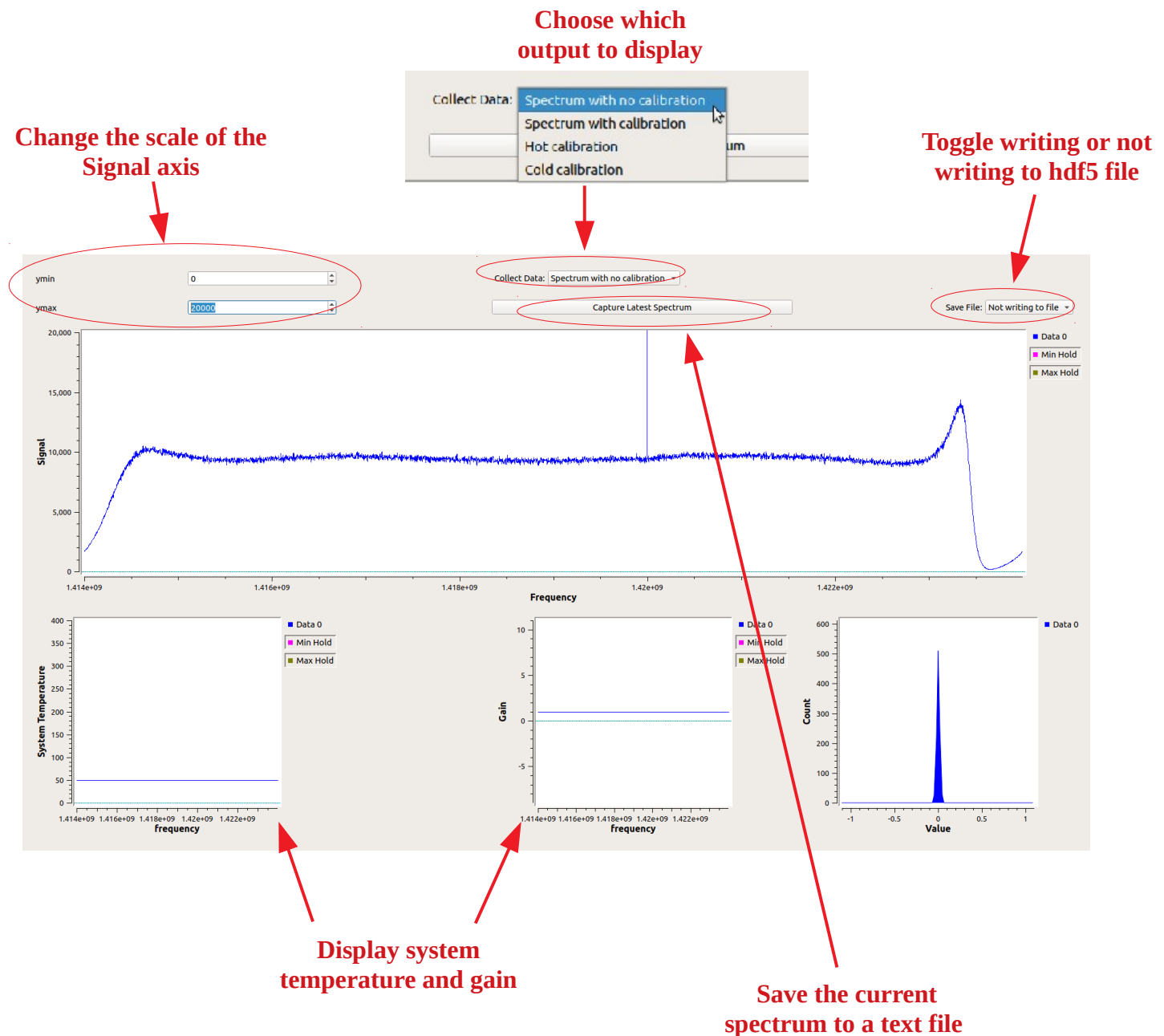

## PROCEDURE: Collecting Spectra Using the *spectrum\_w\_cal.grc* Program

1. Connect horn, LNA, Airspy, and computer USB. Then open the *spectrum\_w\_cal.grc* program in Gnuradio.
2. Before running the program:
  - a. Check the **integration** Variable block for the desired integration time (in seconds). This should be in the upper righthand corner of the Gnuradio screen.
  - b. Check the appropriate file writing Variable blocks, **prefix** and **recfile**. The **prefix** block sets the path to the directory that all files will be written to. The **recfile** is the name of the hdf5 file that will be written.

**NOTE:** Gnuradio does not like folder names with capitals in them. In naming folders, use only lower case, and no spaces in folder names nor in file names!

Other points:

- The spectrum will be output to an hdf5 file (appended by “.h5”) when “**Write to File**” is selected in the **Save File** chooser.

Before running the program, the “**pointing**” and “**notes**” boxes in the hdf5 sink should be changed if you plan on using the hdf5 file.

- i. Whenever the “**Capture Latest Spectrum**” button is pressed, the spectrum in the main spectrum graph will be written to a file in the format: “[timestamp]\_spectrum.csv”. These file names always end in “\_spectrum.csv”, and the recfile name has no bearing on these file names. The timestamp will distinguish between these text files.

The format of these .csv files is *Frequency* in the first column, and *Signal* in the second column. If the spectrum is calibrated, then the Signal is Temperature in Kelvin.

3. Start the program running.
4. Checking the file path and file name: It is advised to run the program for a short time (at least one integration time) and select “Write to File” and hit the “Capture Latest Spectrum” button. Then stop the program, and go to the folder defined in the Prefix box to check that the files were created properly. The misnaming of a file directory path or filename are common errors that will prevent the program from running.
5. When the file paths and names are set, re-start the program.
6. If everything is running okay, the **Count vs. Value** histogram graph in the lower right corner of the display should show some activity by continuously fluctuating. Whenever the program is running, this histogram should continuously fluctuate. If this graph looks frozen, then the program has probably stopped running.
7. Assuming the telescope is pointed toward the sky, a spectrum will appear on the main **Signal vs. Frequency** graph. The first iteration of this graph may take longer than you might expect, up to a minute sometimes, depending on the integration time. Be patient.

The **Collect Data** chooser at the top will be on the **Spectrum with no calibration** option.

The **System Temperature** and **Gain** graphs will display default constant values, as the system has not yet been calibrated.

8. **Calibration:** You are ready to calibrate. Here is the procedure:

- a. Hot calibration:
  - i. Point the horn at the ground at a roughly 45° angle.
  - ii. Change the **Collect Data** chooser at the top to **Hot calibration**.
  - iii. Wait for the graph to settle to a repeated pattern. You may want to adjust the **y<sub>max</sub>** value for the spectrum to fit in the graph.
  - iv. Once you are satisfied that the hot calibration spectrum is steady, switch the **Collect Data** chooser to **Spectrum with no calibration** *before moving the telescope*.
- b. Cold calibration:
  - i. Point the horn at the blank sky where there is no HI signal.
  - ii. Change the **Collect Data** chooser at the top to **Cold calibration**.
  - iii. Wait for the graph to settle to a repeated pattern. You may want to adjust the **y<sub>max</sub>** value.
  - iv. Once you are satisfied that the cold calibration spectrum is steady, switch the **Collect Data** chooser to **Spectrum with no calibration** *before moving the telescope*.
  - v. The **System Temperature** and **Gain** graphs should show some meaningful data as functions of frequency. You are now ready to view calibrated spectra.

## 9. Observations

- a. Select **Spectrum with calibration**. The signal should drop significantly to a value that is less than 100 K. So you will need to change the **y<sub>max</sub>** value accordingly.  
The main graph should be almost completely flat except for where there are HI peaks are.
- b. Depending on how well the system has warmed up, it is possible that the spectrum may drift below 0 K. The data is okay, just offset a little. The **y<sub>min</sub>** value can be changed to negative values to fit the spectrum on the graph if needed.
- c. The active spectrum will be written to the hdf5 file whenever “**Write to File**” is selected in the **Save File** chooser. Data will stop writing to this file when “**Not Writing to File**” is selected, but data will resume being added to the same file if “**Write to File**” is re-selected.
- d. If you wish to update the system temperature and gain, the hot and cold calibrations can be re-done at any time without stopping the program.
- e. Any time you want to write the current spectrum to a “[timestamp]\_\_spectrum.csv” file, hit the “**Capture Latest Spectrum**” button. The timestamp will distinguish these files from one another. It is recommended that you rename them later to more descriptive names.

It is recommended that you keep a log of when the spectra were captured for reference later on when renaming the files.

The format of these .csv files is Frequency in the first column, and Signal in the second column. If the spectrum is calibrated, then the Signal is Temperature in Kelvin.

## 10. **Finish Observations**

- a. When you are ready to quit a run, be sure to select “**Not Writing to File**”, and close the program by hitting “x” in the upper righthand corner. If not, the hdf5 file will not be properly saved.
- b. Check your files and rename accordingly.
- c. Re-starting the program will write to a new hdf5 file with a different time stamp.

## 11. **Recommendations and Tips**

- Plan your viewing ahead of time to help keep track of the .csv files you plan on capturing.
- Keep notes of the files captured to make it easier to sort them out later.
- Capturing the hot and cold calibrations to .csv files can be useful for some post-observation applications.
- CAUTION: If you run the program on a laptop without the power cord plugged in, the computer may go into a sleep mode and stop running the Gnuradio program, depending on the length of the run. Change the power settings on your computer so that it never suspends or goes to sleep.